




**ONE** receiver unit.  
**TWO** tip styles.


**FITS 9 OUT OF 10 PATIENTS**

*WHAT IS easy-CLICK?*


Easy-CLICK is a high definition speaker and ear tip combination that provides exceptional sound quality and high frequency performance.

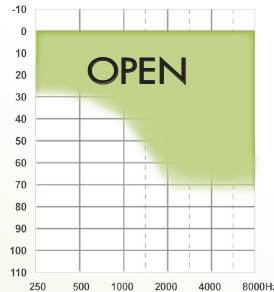

Sizes:  
4mm  
8mm  
10mm

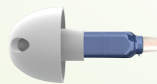


For mild to moderate hearing loss.  
Eliminates the occlusion effect by open style venting.

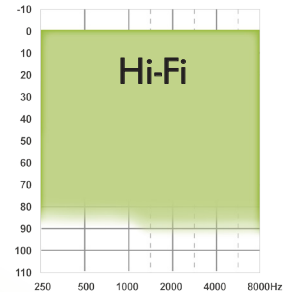

Sizes:  
7mm  
8mm  
9mm  
10mm  
11mm



For mild, moderate or severe hearing loss.  
Eliminates the occlusion effect by using an ultra thin membrane to seal the canal near the second bend.

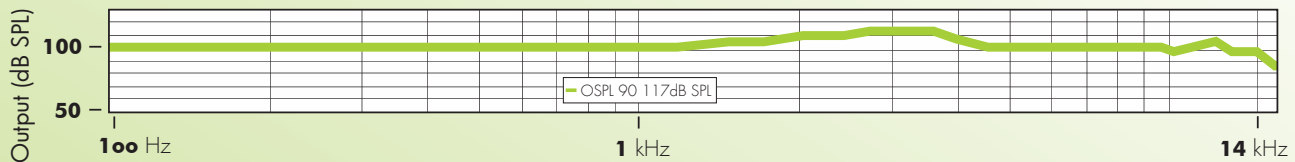


Wideband frequency response from Sebo HD/easy-CLICK hearing instrument

**EASY TO REPLACE**

*DUAL BARRIER WAX FILTER SYSTEM*

Easy-CLICK incorporates a wax filter built into the tip and a second, MPS nano-coated filter on the speaker housing, giving the user unparalleled protection against both wax and moisture. Not only is the new easy-CLICK system more sanitary, it lasts longer and ensures more consistent reliable performance. User-replaceable tips cut down on service visits and repairs saving you both time and money.

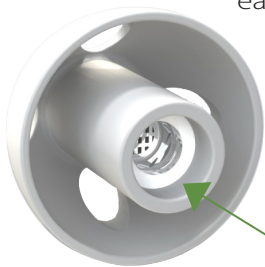

**BUILT IN WAX FILTER**

**AMAZING COMFORT**

*ULTRA SOFT, THIN AND FLEXIBLE*

Easy-CLICK tips utilize an ultra thin membrane for the ultimate comfort experience. Multiple flex zones with 360° mobility ensures the easy-CLICK tip will comfortably adapt vertically and horizontally to canal bends. The ultra thin membrane is made of a memory forming material that retains the true shape of the ear canal after the first use, creating an amazing comfort experience.

